

# ● Analysis Report

## Bruker IVDr Quantification in URine B.I.Quant-UR b<sup>TM</sup>

Sample ID: ALZ\_Urine\_Rack01\_RCM\_221214\_expno800.100000.10r

Measuring Date: 24-Dec-2014 08:45:48

Reporting Date: 12-Dec-2020 09:30:12, 7 page(s), Version 1.1.0

Quantification Method Version: Quant-UR B.1.1.0

### Disclaimer

RESEARCH USE ONLY: This is no clinical diagnostic analysis report. Must not be used for clinical (medical or IVD) diagnosis or for patient management! Additional concentration range information (95% range) provided numerically or graphically in this report must not be used for clinical diagnostic interpretation.

Application of B.I.Quant-UR B 1.1.0 requires use of Bruker's B.I.Methods SOP for urine.

### Summary

The following metabolites were found with concentrations outside the 95% range of Bruker Quant-UR B.1.1.0 urine metabolite concentration database:

Carboxylic acids: *Formic acid (54 mmol/mol Crea), Trigonelline (110 mmol/mol Crea),*


Keto acids and derivatives: *Acetoacetic acid (35 mmol/mol Crea), Oxaloacetic acid (72 mmol/mol Crea).*

Further detailed information is provided on the following pages.

## Contents



|  |          |
|--|----------|
| <b>1 Creatinine</b>                                  | <b>3</b> |
| <b>2 Amines and derivatives</b>                      | <b>3</b> |
| <b>3 Amino acids and derivatives</b>                 | <b>3</b> |
| <b>4 Benzene and substituted derivatives</b>         | <b>4</b> |
| <b>5 Carboxylic acids</b>                            | <b>4</b> |
| <b>6 Fatty acids and derivatives</b>                 | <b>4</b> |
| <b>7 Keto acids and derivatives</b>                  | <b>5</b> |
| <b>8 Purine, Pyridine and Pyrimidine derivatives</b> | <b>5</b> |
| <b>9 Sugars and derivatives</b>                      | <b>5</b> |
| <b>10 Explanations</b>                               | <b>6</b> |

## 1 Creatinine

| Compound   | Conc.<br>mmol/L | LOD<br>mmol/L | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range <sup>(*)</sup><br>mmol/L   |
|------------|-----------------|---------------|-------------|-------------|--------------------|--|
| Creatinine | 3.1             | 0.3           | 3.086       | 100 ●       | 0.074              | 1 - 19  |













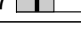
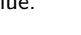
(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 2 Amines and derivatives

| Compound       | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range <sup>(*)</sup><br>$\frac{\text{mmol}}{\text{mol Crea}}$                             |
|----------------|-----------------|--|--|-------------|-------------|--------------------|---|
| Dimethylamine  | 0.14            | 46   | 31   | 0.141       | 100 ●       | 0.006              | $\leq 54$  |
| Trimethylamine | < 0.01          | < 2  | 2  | 0.003       | 96 ●        | 0.001              | $\leq 3$   |




(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 3 Amino acids and derivatives

| Compound             | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range <sup>(*)</sup><br>$\frac{\text{mmol}}{\text{mol Crea}}$                                |
|----------------------|-----------------|--|--|-------------|-------------|--------------------|--|
| 1-Methylhistidine    | < 0.05          | < 15   | 15   | 0.017       | 0 ○         | 0.031              | $\leq 15$   |
| 2-Furoylglycine      | < 0.12          | < 39   | 39   | 0.109       | 99 ●        | 0.010              | $\leq 40$   |
| 4-Aminobutyric acid  | < 0.06          | < 20   | 20   | 0.000       | 0 ○         | 0.450              | $\leq 20$   |
| Alanine              | 0.14            | 45   | 10   | 0.140       | 100 ●       | 0.012              | 11 - 72     |
| Arginine             | < 2.3           | < 750  | 750  | 0.270       | 0 ○         | 1.728              | $\leq 750$  |
| Betaine              | 0.23            | 74   | 7  | 0.227       | 100 ●       | 0.016              | 9 - 78      |
| Creatine             | < 0.15          | < 50   | 50   | 0.088       | 100 ●       | 0.074              | $\leq 280$  |
| Glycine              | 1.2             | 390  | 34   | 1.190       | 100 ●       | 0.032              | 38 - 440    |
| Guanidinoacetic acid | < 0.32          | < 100  | 100  | 0.281       | 98 ●        | 0.050              | $\leq 140$  |
| Methionine           | < 0.06          | < 18   | 18   | 0.000       | 0 ○         | 0.284              | $\leq 18$   |
| N,N-Dimethylglycine  | 0.04            | 14   | 5  | 0.045       | 98 ●        | 0.008              | $\leq 15$   |
| Sarcosine            | 0.01            | 3  | 2  | 0.009       | 41 ○        | 0.006              | $\leq 7$    |
| Taurine              | < 0.44          | < 140  | 140  | 0.204       | 37 ○        | 0.288              | $\leq 170$  |
| Valine               | 0.01            | 4  | 2  | 0.012       | 85 ●        | 0.015              | $\leq 7$    |









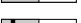

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 4 Benzene and substituted derivatives

| Compound        | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range(*)<br>$\frac{\text{mmol}}{\text{mol Crea}}$                                      |
|-----------------|-----------------|--|--|-------------|-------------|--------------------|--|
| Benzoic acid    | < 0.03          | < 10   | 10   | 0.000       | 0○          | 0.017              | ≤ 10    |
| D-Mandelic acid | < 0.01          | < 2  | 2  | 0.000       | 0○          | 0.006              | 2 - 17  |
| Hippuric acid   | 1.2             | 370  | 170  | 1.152       | 100●        | 0.062              | ≤ 660   |


(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 5 Carboxylic acids

| Compound        | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range(*)<br>$\frac{\text{mmol}}{\text{mol Crea}}$                                       |
|-----------------|-----------------|--|--|-------------|-------------|--------------------|---|
| Acetic acid     | 0.05            | 16   | 5  | 0.049       | 78○         | 0.036              | ≤ 51     |
| Citric acid     | 1.8             | 580  | 40   | 1.777       | 100●        | 0.171              | ≤ 700   |
| Formic acid     | 0.17            | 54   | 10   | 0.166       | 100●        | 0.008              | ≤ 43   |
| Fumaric acid    | < 0.01          | < 2  | 2  | 0.004       | 99●         | 0.000              | ≤ 3    |
| Imidazole       | < 0.15          | < 48   | 48   | 0.000       | 0○          | 0.110              | ≤ 48   |
| Lactic acid     | < 0.15          | < 49   | 49   | 0.059       | 93●         | 0.094              | ≤ 110  |
| Proline betaine | < 0.08          | < 25   | 25   | 0.072       | 71○         | 0.068              | ≤ 280  |
| Succinic acid   | 0.04            | 12   | 5  | 0.036       | 69○         | 0.019              | ≤ 39   |
| Tartaric acid   | < 0.01          | < 5  | 5  | 0.007       | 41○         | 0.010              | ≤ 110  |
| Trigonelline    | 0.34            | 110  | 35   | 0.337       | 100●        | 0.007              | ≤ 67   |






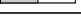
(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 6 Fatty acids and derivatives

| Compound              | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range(*)<br>$\frac{\text{mmol}}{\text{mol Crea}}$                                      |
|-----------------------|-----------------|--|--|-------------|-------------|--------------------|--|
| 2-Methylsuccinic acid | < 0.15          | < 48   | 48   | 0.000       | 0○          | 0.207              | ≤ 48  |






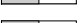

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 7 Keto acids and derivatives

| Compound              | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range(*)<br>$\frac{\text{mmol}}{\text{mol Crea}}$                                     |
|-----------------------|-----------------|--|--|-------------|-------------|--------------------|---|
| 2-Oxoglutaric acid    | < 0.28          | < 92   | 92   | 0.086       | 57 ○        | 0.190              | ≤ 92   |
| 3-Hydroxybutyric acid | < 0.32          | < 100  | 100  | 0.000       | 0 ○         | 0.883              | ≤ 100  |
| Acetoacetic acid      | 0.11            | 35   | 14   | 0.108       | 98 ●        | 0.028              | ≤ 30   |
| Acetone               | 0.02            | 5  | 2  | 0.016       | 97 ●        | 0.002              | ≤ 7    |
| Oxaloacetic acid      | 0.22            | 72   | 17   | 0.223       | 94 ●        | 0.082              | ≤ 66   |
| Pyruvic acid          | < 0.03          | < 9  | 9  | 0.016       | 89 ●        | 0.005              | ≤ 13   |





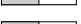

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 8 Purine, Pyridine and Pyrimidine derivatives

| Compound             | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range(*)<br>$\frac{\text{mmol}}{\text{mol Crea}}$                                       |
|----------------------|-----------------|--|--|-------------|-------------|--------------------|---|
| 1-Methyladenosine    | < 0.01          | < 5  | 5  | 0.000       | 0 ○         | 0.087              | ≤ 5    |
| 1-Methylnicotinamide | < 0.10          | < 32   | 32   | 0.026       | 90 ●        | 0.011              | ≤ 32   |
| Adenosine            | < 1.2           | < 390  | 390  | 0.000       | 0 ○         | 0.879              | ≤ 390  |
| Allantoin            | < 0.05          | < 17   | 17   | 0.023       | 80 ○        | 0.011              | ≤ 47   |
| Allopurinol          | < 0.03          | < 10   | 10   | 0.016       | 74 ○        | 0.023              | ≤ 11   |
| Caffeine             | < 0.14          | < 45   | 45   | 0.102       | 66 ○        | 0.134              | ≤ 61   |
| Inosine              | < 0.06          | < 19   | 19   | 0.021       | 72 ○        | 0.027              | ≤ 19   |

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 9 Sugars and derivatives

| Compound     | Conc.<br>mmol/L | Conc.<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | LOD<br>$\frac{\text{mmol}}{\text{mol Crea}}$ | r<br>mmol/L | $\rho$<br>% | $\Delta$<br>mmol/L | 95% Range(*)<br>$\frac{\text{mmol}}{\text{mol Crea}}$  |
|--------------|-----------------|--|--|-------------|-------------|--------------------|--|
| D-Galactose  | < 0.13          | < 43   | 43   | 0.000       | 0 ○         | 0.012              | ≤ 44    |
| D-Glucose    | 0.24            | 77   | 34   | 0.238       | 87 ●        | 0.053              | ≤ 140   |
| D-Lactose    | < 0.30          | < 96   | 96   | 0.090       | 89 ●        | 0.027              | ≤ 96    |
| D-Mannitol   | < 0.57          | < 180  | 180  | 0.000       | 0 ○         | 1.596              | ≤ 180   |
| D-Mannose    | < 0.02          | < 6  | 6  | 0.000       | 0 ○         | 0.054              | ≤ 8     |
| Myo-Inositol | < 14            | < 4400   | 4400   | 0.000       | 0 ○         | 3.359              | ≤ 4400  |

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 10 Explanations

This section contains the definition of the parameters used above. In the section 10.1 a short manual, how to interpret the results, is presented. The section 10.3 contains the exact definitions of the parameters  $r$ ,  $\rho$  and  $\Delta$ .

### 10.1 How to read the result



Figure 1: Examples of fitting.

In the figure 1(a), the black line, the blue line and the yellow line represent the original spectrum, the calculated signal fit and its baseline, respectively.

The blue area relates to the metabolite concentration to be determined and the red area represents a residue.

In case of the signal overlap a different approach is used: two or more overlapping signals are being fitted simultaneously. The most iconic example of such signals are the ones generated by  $\text{CH}_3$  groups of Creatinine and Creatine. In such a case, the blue line and the grey area relate the sum of all fitted signals. The blue area corresponds to the concentration of the metabolite of interest (cf. figure 1(b)).

### 10.2 Result parameters

- Conc.** is the final result concentration of the metabolite,
- LOD** is the *limit of detection* of the given metabolite,
- r** is the *raw concentration* i.e. the concentration equivalent of the resulting signal fit prior to comparing to **LOD** (relates to the blue area, cf.  $\alpha$ ),
- $\rho$  is the correlation of lineshape metabolite signal with calculated fit characterizing the match between metabolite signal and fit (cf.  $\beta$ ). Depending on the value of  $\rho$ , the following *flag* is displayed:

- ●, if the correlation is 95%,
  - ●, if the correlation is in between 85% and 95%,
  - ○, if the correlation is less than 85%,
- e)  $\Delta$  is the concentration equivalent of the difference between metabolite signal and calculated fit (residue corresponding to the **the red area**, cf.  $\gamma$ )).

### 10.3 Detailed definitions

Let  $s$ ,  $f$  and  $b$  denote the functions describing the *raw spectra*, *fitted curve* and *(fitted) baseline* respectively. These functions are chosen such that  $s \approx f + b$ . Moreover, let  $I$  be a relevant PPM interval and  $P_N$  be the proton number for given metabolite/signal.

$\alpha$ )  $r$  (*raw concentration*) is defined as

$$r = \frac{1}{P_N} \int_{\mathbb{R}} f(\xi) d\xi.$$

$\beta$ )  $\rho$  is the *correlation* of the functions  $s$  and  $f + b$ , i.e.

$$\rho = \max(0, \text{corr}(\bar{s}, \overline{f+b})),$$

where  $\bar{s}$ ,  $\overline{f+b}$  are numerical representations of the functions  $s$  and  $f + b$  on sufficiently fine mesh of the interval  $I$ .

$\gamma$ )  $\Delta$  is the the area between the raw signal  $s$  and the fitted data  $f + b$  on the interval  $I$  expressed in the terms of the concentration, i.e.

$$\Delta = \frac{1}{P_N} \int_I |s(\xi) - f(\xi) - b(\xi)| d\xi.$$