

# ● Analysis Report

## Bruker IVDr Quantification in URine B.I.Quant-UR b<sup>TM</sup>

Sample ID: ALZ\_Urine\_Rack01\_RCM\_221214\_expno550.100000.10r

Measuring Date: 24-Dec-2014 02:59:48

Reporting Date: 11-Dec-2020 21:53:03, 7 page(s), Version 1.1.0

Quantification Method Version: Quant-UR B.1.1.0

### Disclaimer

RESEARCH USE ONLY: This is no clinical diagnostic analysis report. Must not be used for clinical (medical or IVD) diagnosis or for patient management! Additional concentration range information (95% range) provided numerically or graphically in this report must not be used for clinical diagnostic interpretation.

Application of B.I.Quant-UR B 1.1.0 requires use of Bruker's B.I.Methods SOP for urine.

### Summary

The following metabolites were found with concentrations outside the 95% range of Bruker Quant-UR B.1.1.0 urine metabolite concentration database:


*Benzene and substituted derivatives: Benzoic acid (47 mmol/mol Crea).*

Further detailed information is provided on the following pages.

## Contents



<b>1 Creatinine</b>	<b>3</b>
<b>2 Amines and derivatives</b>	<b>3</b>
<b>3 Amino acids and derivatives</b>	<b>3</b>
<b>4 Benzene and substituted derivatives</b>	<b>4</b>
<b>5 Carboxylic acids</b>	<b>4</b>
<b>6 Fatty acids and derivatives</b>	<b>4</b>
<b>7 Keto acids and derivatives</b>	<b>5</b>
<b>8 Purine, Pyridine and Pyrimidine derivatives</b>	<b>5</b>
<b>9 Sugars and derivatives</b>	<b>5</b>
<b>10 Explanations</b>	<b>6</b>

## 1 Creatinine

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range <sup>(*)</sup> mmol/L
Creatinine	3.9	0.3	3.923	100 ●	0.052	1 - 19 













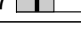
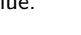
(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 2 Amines and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range <sup>(*)</sup> $\frac{\text{mmol}}{\text{mol Crea}}$
Dimethylamine	0.14	36	31	0.142	100 ●	0.003	$\leq 54$ 
Trimethylamine	< 0.01	< 2	2	0.003	92 ●	0.001	$\leq 3$ 

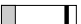


(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 3 Amino acids and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range <sup>(*)</sup> $\frac{\text{mmol}}{\text{mol Crea}}$
1-Methylhistidine	< 0.06	< 15	15	0.000	0 ○	0.044	$\leq 15$ 
2-Furoylglycine	< 0.15	< 39	39	0.055	96 ●	0.008	$\leq 40$ 
4-Aminobutyric acid	< 0.08	< 20	20	0.000	0 ○	0.152	$\leq 20$ 
Alanine	0.13	32	10	0.127	100 ●	0.005	11 - 72 
Arginine	< 2.9	< 750	750	0.257	19 ○	1.791	$\leq 750$ 
Betaine	0.05	12	7	0.047	100 ●	0.009	9 - 78 
Creatine	< 0.20	< 50	50	0.074	100 ●	0.052	$\leq 280$ 
Glycine	0.54	140	34	0.539	100 ●	0.015	38 - 440 
Guanidinoacetic acid	< 0.41	< 100	100	0.201	94 ●	0.052	$\leq 140$ 
Methionine	< 0.07	< 18	18	0.000	0 ○	0.493	$\leq 18$ 
N,N-Dimethylglycine	0.04	10	5	0.041	97 ●	0.007	$\leq 15$ 
Sarcosine	< 0.01	< 2	2	0.003	0 ○	0.007	$\leq 7$ 
Taurine	0.60	150	140	0.604	94 ●	0.223	$\leq 170$ 
Valine	0.01	4	2	0.015	5 ○	0.021	$\leq 7$ 





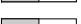





(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 4 Benzene and substituted derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
Benzoic acid	0.19	47	10	0.186	99 ●	0.018	$\leq 10$ 
D-Mandelic acid	< 0.01	< 2	2	0.000	0 ○	0.071	2 - 17 
Hippuric acid	< 0.67	< 170	170	0.582	98 ●	0.100	$\leq 660$ 


(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 5 Carboxylic acids

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
Acetic acid	0.08	20	5	0.079	99 ●	0.009	$\leq 51$ 
Citric acid	1.6	400	40	1.557	100 ●	0.152	$\leq 700$ 
Formic acid	0.09	22	10	0.086	100 ●	0.004	$\leq 43$ 
Fumaric acid	0.01	2	2	0.009	99 ●	0.001	$\leq 3$ 
Imidazole	< 0.19	< 48	48	0.000	0 ○	0.138	$\leq 48$ 
Lactic acid	< 0.19	< 49	49	0.126	96 ●	0.029	$\leq 110$ 
Proline betaine	< 0.10	< 25	25	0.040	62 ○	0.068	$\leq 280$ 
Succinic acid	0.07	17	5	0.067	97 ●	0.012	$\leq 39$ 
Tartaric acid	< 0.02	< 5	5	0.012	29 ○	0.018	$\leq 110$ 
Trigonelline	0.25	65	35	0.253	100 ●	0.006	$\leq 67$ 






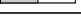
(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 6 Fatty acids and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
2-Methylsuccinic acid	< 0.19	< 48	48	0.033	0 ○	0.045	$\leq 48$ 






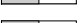

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 7 Keto acids and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
2-Oxoglutaric acid	< 0.36	< 92	92	0.110	82 ○	0.172	≤ 92 
3-Hydroxybutyric acid	< 0.40	< 100	100	0.000	0 ○	0.363	≤ 100 
Acetoacetic acid	< 0.06	< 14	14	0.051	87 ●	0.025	≤ 30 
Acetone	0.01	3	2	0.013	97 ●	0.003	≤ 7 
Oxaloacetic acid	0.07	17	17	0.067	84 ○	0.055	≤ 66 
Pyruvic acid	< 0.04	< 9	9	0.015	89 ●	0.006	≤ 13 





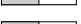

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 8 Purine, Pyridine and Pyrimidine derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
1-Methyladenosine	< 0.02	< 5	5	0.000	0 ○	0.089	≤ 5 
1-Methylnicotinamide	< 0.12	< 32	32	0.026	95 ●	0.008	≤ 32 
Adenosine	< 1.5	< 390	390	0.000	0 ○	0.964	≤ 390 
Allantoin	< 0.07	< 17	17	0.043	99 ●	0.003	≤ 47 
Allopurinol	< 0.04	< 10	10	0.027	89 ●	0.042	≤ 11 
Caffeine	< 0.18	< 45	45	0.095	83 ○	0.087	≤ 61 
Inosine	< 0.07	< 19	19	0.064	51 ○	0.061	≤ 19 

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 9 Sugars and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
D-Galactose	< 0.17	< 43	43	0.013	97 ●	0.001	≤ 44 
D-Glucose	< 0.13	< 34	34	0.130	63 ○	0.126	≤ 140 
D-Lactose	< 0.38	< 96	96	0.100	93 ●	0.024	≤ 96 
D-Mannitol	< 0.72	< 180	180	0.000	0 ○	1.527	≤ 180 
D-Mannose	< 0.02	< 6	6	0.000	0 ○	0.057	≤ 8 
Myo-Inositol	< 17	< 4400	4400	0.000	0 ○	2.965	≤ 4400 

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 10 Explanations

This section contains the definition of the parameters used above. In the section 10.1 a short manual, how to interpret the results, is presented. The section 10.3 contains the exact definitions of the parameters  $r$ ,  $\rho$  and  $\Delta$ .

### 10.1 How to read the result

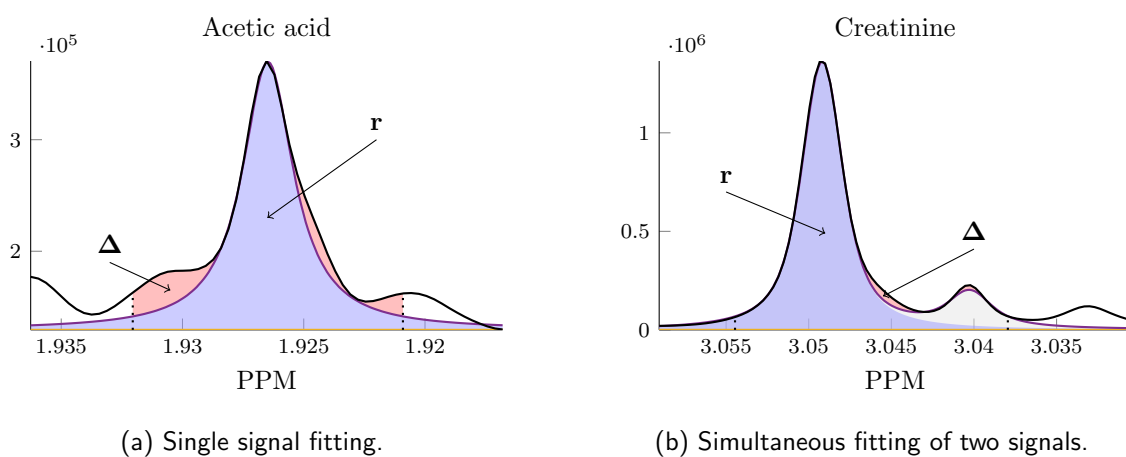


Figure 1: Examples of fitting.

In the figure 1(a), the black line, the blue line and the yellow line represent the original spectrum, the calculated signal fit and its baseline, respectively.

The blue area relates to the metabolite concentration to be determined and the red area represents a residue.

In case of the signal overlap a different approach is used: two or more overlapping signals are being fitted simultaneously. The most iconic example of such signals are the ones generated by  $\text{CH}_3$  groups of Creatinine and Creatine. In such a case, the blue line and the grey area relate the sum of all fitted signals. The blue area corresponds to the concentration of the metabolite of interest (cf. figure 1(b)).

### 10.2 Result parameters

- Conc.** is the final result concentration of the metabolite,
- LOD** is the *limit of detection* of the given metabolite,
- r** is the *raw concentration* i.e. the concentration equivalent of the resulting signal fit prior to comparing to **LOD** (relates to the blue area, cf.  $\alpha$ ),
- $\rho$  is the correlation of lineshape metabolite signal with calculated fit characterizing the match between metabolite signal and fit (cf.  $\beta$ ). Depending on the value of  $\rho$ , the following *flag* is displayed:

- ●, if the correlation is 95%,
  - ●, if the correlation is in between 85% and 95%,
  - ○, if the correlation is less than 85%,
- e)  $\Delta$  is the concentration equivalent of the difference between metabolite signal and calculated fit (residue corresponding to the the red area, cf.  $\gamma$ )).

### 10.3 Detailed definitions

Let  $s$ ,  $f$  and  $b$  denote the functions describing the *raw spectra*, *fitted curve* and *(fitted) baseline* respectively. These functions are chosen such that  $s \approx f + b$ . Moreover, let  $I$  be a relevant PPM interval and  $P_N$  be the proton number for given metabolite/signal.

$\alpha$ )  $r$  (*raw concentration*) is defined as

$$r = \frac{1}{P_N} \int_{\mathbb{R}} f(\xi) d\xi.$$

$\beta$ )  $\rho$  is the *correlation* of the functions  $s$  and  $f + b$ , i.e.

$$\rho = \max(0, \text{corr}(\bar{s}, \overline{f+b})),$$

where  $\bar{s}$ ,  $\overline{f+b}$  are numerical representations of the functions  $s$  and  $f + b$  on sufficiently fine mesh of the interval  $I$ .

$\gamma$ )  $\Delta$  is the the area between the raw signal  $s$  and the fitted data  $f + b$  on the interval  $I$  expressed in the terms of the concentration, i.e.

$$\Delta = \frac{1}{P_N} \int_I |s(\xi) - f(\xi) - b(\xi)| d\xi.$$