

# ● Analysis Report

## Bruker IVDr Quantification in URine B.I.Quant-UR b<sup>TM</sup>

Sample ID: ALZ\_Urine\_Rack01\_RCM\_221214\_expno320.100000.10r

Measuring Date: 23-Dec-2014 21:21:19

Reporting Date: 12-Dec-2020 15:24:22, 7 page(s), Version 1.1.0

Quantification Method Version: Quant-UR B.1.1.0

### Disclaimer

RESEARCH USE ONLY: This is no clinical diagnostic analysis report. Must not be used for clinical (medical or IVD) diagnosis or for patient management! Additional concentration range information (95% range) provided numerically or graphically in this report must not be used for clinical diagnostic interpretation.

Application of B.I.Quant-UR B 1.1.0 requires use of Bruker's B.I.Methods SOP for urine.

### Summary

The following metabolites were found with concentrations outside the 95% range of Bruker Quant-UR B.1.1.0 urine metabolite concentration database:


Amino acids and derivatives: 4-Aminobutyric acid (21 mmol/mol Crea).

Further detailed information is provided on the following pages.

## Contents


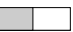
<b>1 Creatinine</b>	<b>3</b>
<b>2 Amines and derivatives</b>	<b>3</b>
<b>3 Amino acids and derivatives</b>	<b>3</b>
<b>4 Benzene and substituted derivatives</b>	<b>4</b>
<b>5 Carboxylic acids</b>	<b>4</b>
<b>6 Fatty acids and derivatives</b>	<b>4</b>
<b>7 Keto acids and derivatives</b>	<b>5</b>
<b>8 Purine, Pyridine and Pyrimidine derivatives</b>	<b>5</b>
<b>9 Sugars and derivatives</b>	<b>5</b>
<b>10 Explanations</b>	<b>6</b>

## 1 Creatinine

Compound	Conc. mmol/L	LOD mmol/L	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range <sup>(*)</sup> mmol/L
Creatinine	8.6	0.3	8.558	100 ●	0.085	1 - 19 















<sup>(\*)</sup> Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 2 Amines and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range <sup>(*)</sup> $\frac{\text{mmol}}{\text{mol Crea}}$
Dimethylamine	0.29	34	31	0.287	100 ●	0.010	$\leq 54$ 
Trimethylamine	< 0.02	< 2	2	0.001	0 ○	0.003	$\leq 3$ 

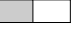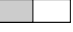


<sup>(\*)</sup> Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 3 Amino acids and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range <sup>(*)</sup> $\frac{\text{mmol}}{\text{mol Crea}}$
1-Methylhistidine	< 0.13	< 15	15	0.000	0 ○	0.217	$\leq 15$ 
2-Furoylglycine	< 0.33	< 39	39	0.111	49 ○	0.039	$\leq 40$ 
4-Aminobutyric acid	0.18	21	20	0.180	48 ○	0.212	$\leq 20$ 
Alanine	0.27	32	10	0.275	100 ●	0.022	11 - 72 
Arginine	< 6.4	< 750	750	0.694	0 ○	2.658	$\leq 750$ 
Betaine	0.20	23	7	0.197	100 ●	0.030	9 - 78 
Creatine	< 0.43	< 50	50	0.275	100 ●	0.085	$\leq 280$ 
Glycine	0.93	110	34	0.934	100 ●	0.038	38 - 440 
Guanidinoacetic acid	< 0.88	< 100	100	0.356	92 ●	0.110	$\leq 140$ 
Methionine	< 0.15	< 18	18	0.000	0 ○	0.445	$\leq 18$ 
N,N-Dimethylglycine	< 0.04	< 5	5	0.043	98 ●	0.005	$\leq 15$ 
Sarcosine	< 0.02	< 2	2	0.000	0 ○	0.005	$\leq 7$ 
Taurine	< 1.2	< 140	140	0.720	86 ●	0.278	$\leq 170$ 
Valine	0.03	4	2	0.033	95 ●	0.007	$\leq 7$ 











<sup>(\*)</sup> Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 4 Benzene and substituted derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
Benzoic acid	< 0.08	< 10	10	0.027	93 ●	0.007	≤ 10 
D-Mandelic acid	< 0.02	< 2	2	0.000	0 ○	0.111	2 - 17 
Hippuric acid	3.5	400	170	3.458	99 ●	0.351	≤ 660 


(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 5 Carboxylic acids

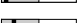Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
Acetic acid	0.20	24	5	0.203	99 ●	0.020	≤ 51 
Citric acid	2.4	280	40	2.361	100 ●	0.289	≤ 700 
Formic acid	0.28	32	10	0.276	100 ●	0.017	≤ 43 
Fumaric acid	< 0.02	< 2	2	0.005	99 ●	0.000	≤ 3 
Imidazole	< 0.41	< 48	48	0.078	3 ○	0.077	≤ 48 
Lactic acid	< 0.42	< 49	49	0.228	99 ●	0.033	≤ 110 
Proline betaine	0.30	35	25	0.304	99 ●	0.028	≤ 280 
Succinic acid	0.11	13	5	0.107	98 ●	0.014	≤ 39 
Tartaric acid	0.13	16	5	0.135	99 ●	0.012	≤ 110 
Trigonelline	0.34	39	35	0.337	100 ●	0.009	≤ 67 






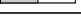
(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 6 Fatty acids and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
2-Methylsuccinic acid	< 0.41	< 48	48	0.051	0 ○	0.110	≤ 48 






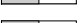

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 7 Keto acids and derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
2-Oxoglutaric acid	< 0.79	< 92	92	0.013	0 ○	0.153	≤ 92 
3-Hydroxybutyric acid	< 0.88	< 100	100	0.102	26 ○	0.357	≤ 100 
Acetoacetic acid	0.12	14	14	0.120	94 ●	0.053	≤ 30 
Acetone	0.02	2	2	0.020	99 ●	0.005	≤ 7 
Oxaloacetic acid	0.19	22	17	0.191	89 ●	0.138	≤ 66 
Pyruvic acid	< 0.08	< 9	9	0.038	96 ●	0.008	≤ 13 







(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 8 Purine, Pyridine and Pyrimidine derivatives

Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
1-Methyladenosine	< 0.04	< 5	5	0.000	0 ○	0.144	≤ 5 
1-Methylnicotinamide	< 0.27	< 32	32	0.057	99 ●	0.004	≤ 32 
Adenosine	< 3.3	< 390	390	0.000	0 ○	1.696	≤ 390 
Allantoin	< 0.14	< 17	17	0.092	98 ●	0.011	≤ 47 
Allopurinol	< 0.09	< 10	10	0.048	76 ○	0.048	≤ 11 
Caffeine	< 0.39	< 45	45	0.185	87 ●	0.215	≤ 61 
Inosine	< 0.16	< 19	19	0.032	92 ●	0.049	≤ 19 

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 9 Sugars and derivatives

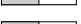Compound	Conc. mmol/L	Conc. $\frac{\text{mmol}}{\text{mol Crea}}$	LOD $\frac{\text{mmol}}{\text{mol Crea}}$	r mmol/L	$\rho$ %	$\Delta$ mmol/L	95% Range(*) $\frac{\text{mmol}}{\text{mol Crea}}$
D-Galactose	< 0.37	< 43	43	0.033	91 ●	0.004	≤ 44 
D-Glucose	0.79	93	34	0.793	93 ●	0.134	≤ 140 
D-Lactose	< 0.82	< 96	96	0.202	85 ●	0.091	≤ 96 
D-Mannitol	< 1.6	< 180	180	0.000	0 ○	3.909	≤ 180 
D-Mannose	< 0.05	< 6	6	0.031	0 ○	0.082	≤ 8 
Myo-Inositol	< 38	< 4400	4400	0.000	0 ○	10.34	≤ 4400 

(\*) Gray horizontal boxes represent 95% concentration range, black vertical lines represent sample value.

## 10 Explanations

This section contains the definition of the parameters used above. In the section 10.1 a short manual, how to interpret the results, is presented. The section 10.3 contains the exact definitions of the parameters  $r$ ,  $\rho$  and  $\Delta$ .

### 10.1 How to read the result

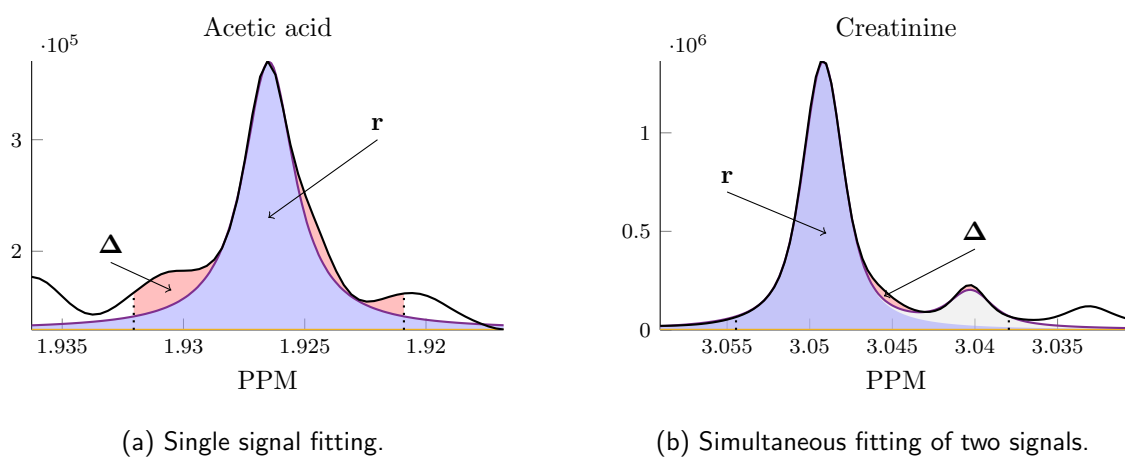


Figure 1: Examples of fitting.

In the figure 1(a), the black line, the blue line and the yellow line represent the original spectrum, the calculated signal fit and its baseline, respectively.

The blue area relates to the metabolite concentration to be determined and the red area represents a residue.

In case of the signal overlap a different approach is used: two or more overlapping signals are being fitted simultaneously. The most iconic example of such signals are the ones generated by  $\text{CH}_3$  groups of Creatinine and Creatine. In such a case, the blue line and the grey area relate the sum of all fitted signals. The blue area corresponds to the concentration of the metabolite of interest (cf. figure 1(b)).

### 10.2 Result parameters

- Conc.** is the final result concentration of the metabolite,
- LOD** is the *limit of detection* of the given metabolite,
- r** is the *raw concentration* i.e. the concentration equivalent of the resulting signal fit prior to comparing to **LOD** (relates to the blue area, cf.  $\alpha$ ),
- $\rho$  is the correlation of lineshape metabolite signal with calculated fit characterizing the match between metabolite signal and fit (cf.  $\beta$ ). Depending on the value of  $\rho$ , the following *flag* is displayed:

- ●, if the correlation is 95%,
  - ●, if the correlation is in between 85% and 95%,
  - ○, if the correlation is less than 85%,
- e)  $\Delta$  is the concentration equivalent of the difference between metabolite signal and calculated fit (residue corresponding to the **the red area**, cf.  $\gamma$ )).

### 10.3 Detailed definitions

Let  $s$ ,  $f$  and  $b$  denote the functions describing the *raw spectra*, *fitted curve* and *(fitted) baseline* respectively. These functions are chosen such that  $s \approx f + b$ . Moreover, let  $I$  be a relevant PPM interval and  $P_N$  be the proton number for given metabolite/signal.

$\alpha$ )  $r$  (*raw concentration*) is defined as

$$r = \frac{1}{P_N} \int_{\mathbb{R}} f(\xi) d\xi.$$

$\beta$ )  $\rho$  is the *correlation* of the functions  $s$  and  $f + b$ , i.e.

$$\rho = \max(0, \text{corr}(\bar{s}, \overline{f+b})) ,$$

where  $\bar{s}$ ,  $\overline{f+b}$  are numerical representations of the functions  $s$  and  $f + b$  on sufficiently fine mesh of the interval  $I$ .

$\gamma$ )  $\Delta$  is the the area between the raw signal  $s$  and the fitted data  $f + b$  on the interval  $I$  expressed in the terms of the concentration, i.e.

$$\Delta = \frac{1}{P_N} \int_I |s(\xi) - f(\xi) - b(\xi)| d\xi.$$